

AeraSeal™

Designed to facilitate air and CO₂ exchange

Application:

Cell, tissue and bacterial cultures


Temperature Range:

-20°C to +80°C

CATALOG #	DESCRIPTION	QUANTITY
B-100	AeraSeal	100
BS-25	AeraSeal, Sterile	50

AeraSeal features a hydrophobic, gas permeable membrane with a medical grade, non-cytotoxic adhesive. The porous membrane consists of layers of heat welded rayon fibers.

- Pore size from 10 µm to 50 µm
- Moisture vapor transmission rate = 4,200g/m²/24hr
- Easy-to-pierce with serological pipets and pipet tips for direct sample recovery
- Available in roll format for automation: AeraSeal Roll-Seal

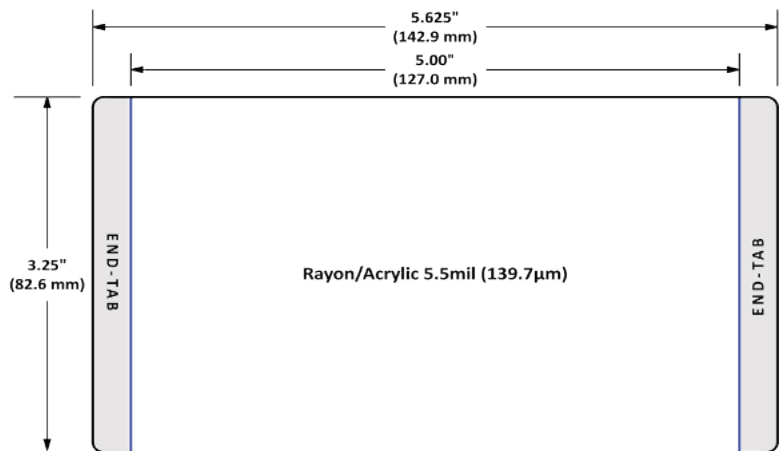

"For a faster and easier way to apply via rolls and a dispenser, AeraSeal System"

