



Zcap™ 1x8 PCR FoilStrips™

For sealing a PCR 8-tube strip or one row of a 96-well plate

Application:

PCR and light sensitive assays


Temperature Range:

-20°C to +120°C

CATALOG #	DESCRIPTION	QUANTITY
F2-CS1X8-300	EZcap 1x8 PCR FoilStrips	300

An economical and more user-friendly alternative to PCR strip caps; no sore thumbs and fingers guaranteed! Precision cut to a smaller dimension for sealing one PCR 8-tube strip or a single row of a 96-well PCR or deep-well storage plate whenever rows must be selectively protected or accessed. FoilStrips remove easily with no splash back as can occur with strip-caps.

- Certified DNase- & RNase- free
- Each sheet contains 6 strips on a continuous liner for ease of use
- Pierceable with a pipet tip or robotic probe for direct sample recovery
- The conformable nature of foil allows for the visual verification of the seal's security around each well
- 2.0mil aluminum foil removes easily without tearing

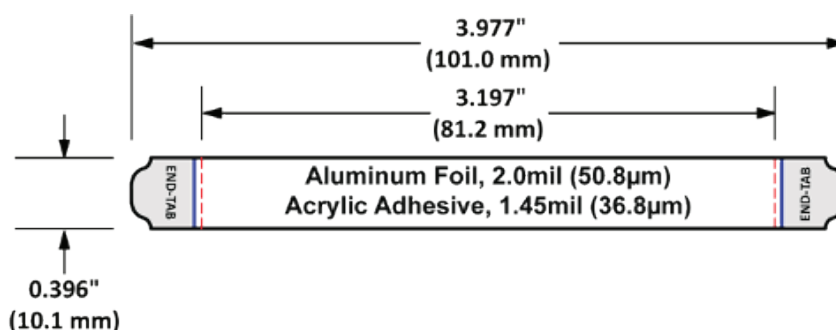

"For a full plate size version of this foil, check out FoilSeal"

