

# SealMate™ AeraSeal™


For fast and efficient microplate sealing

## Application:

Cell, tissue and bacterial cultures

## Temperature Range:



-20°C to +80°C

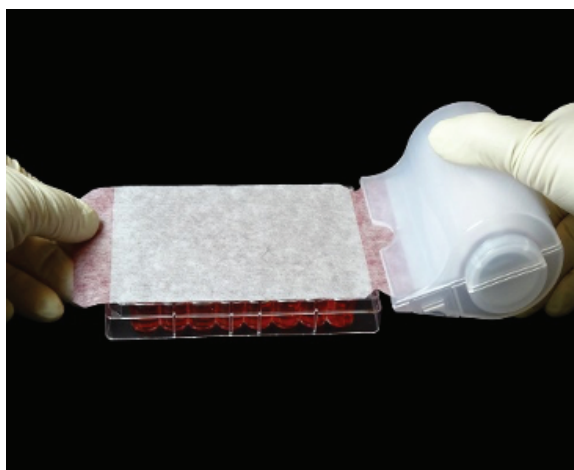

CATALOG #	DESCRIPTION	QUANTITY
SM-KIT-B	AeraSeal Starter Kit, Dispenser + 2 Rolls	50/Roll
SM-KIT-BS	AeraSeal Starter Kit, Dispenser + 2 Rolls, Sterile	50/Roll
RSMB-2	AeraSeal Film, 2 Rolls	50/Roll
RSMB-2-S	AeraSeal Film, 2 Rolls, Sterile	50/Roll

The SealMate AeraSeal Starter Kit features a re-usable dispenser and two rolls of perforated film. The dispenser facilitates the quick and easy dispensing of film from a roll allowing microplates to be sealed three times faster than using sheeted films. Cost-effective replacement rolls are sold separately. Ideal for use under the hood; the SealMate dispenser minimizes film contamination.

- Environmentally friendly; no paper liner to remove and discard; utilizes approximately 59% less material than AeraSeal sheeted film
- Each sheet is separated by an adhesive-free, perforated, color-coded red band
- The smooth lower surface of the dispenser can be used to seal the film to the plate
- Dispenser protects unused seals from contamination
- Lower support prop prevents dispensing bill from contacting surfaces to prevent contamination
- Replacement rolls can be easily installed by releasing the locking clasps on the dispenser
- Designed to facilitate air and CO<sub>2</sub> exchange
- Pore size from 10µm to 50µm
- Moisture vapor transmission rate = 4,200g/m<sup>2</sup>/24hr
- Easy-to-pierce with serological pipets and pipet tips for direct sample recovery

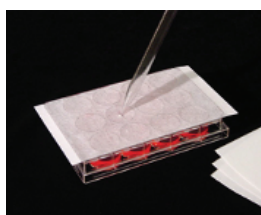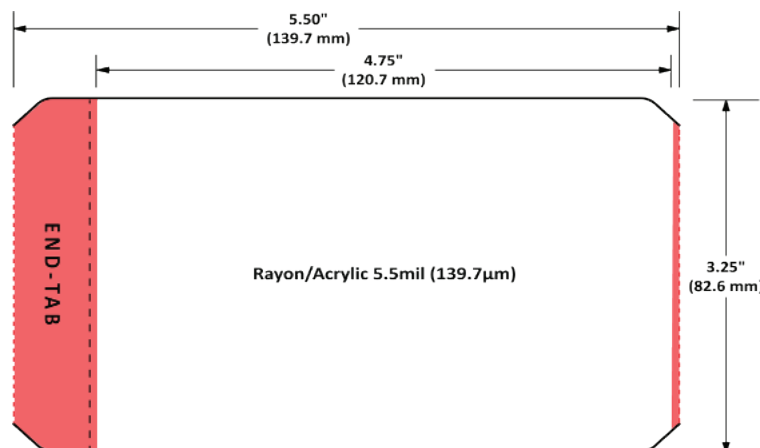

*"For the sheeted alternative,  
Check out AeraSeal individual sheets."*