

Zone-Free™

For easy sample access without adhesive fouling of pipet tips or probes

Application:

Short-term storage and robotics



Temperature Range:

-40°C to +90°C

CATALOG #	DESCRIPTION	QUANTITY
ZAF-PE-50	Zone-Free	50

Zone-Free is a multilayer film consisting of a top pierceable clear polyethylene film and a white polypropylene sublayer with 96 circular areas void of adhesive. The white sublayer aids in the alignment of the film to the plate; alignment of opposing corner wells ensures that all 96 wells will be properly aligned. The adhesive-free zones are 9mm on center to facilitate multi-channel pipetting without adhesive fouling of pipet tips or robotic probes.

- Pierceable
- Offers excellent chemical resistance
- Samples have minimal contact with adhesive due to the adhesive-free zones

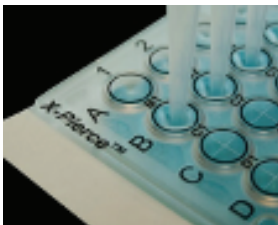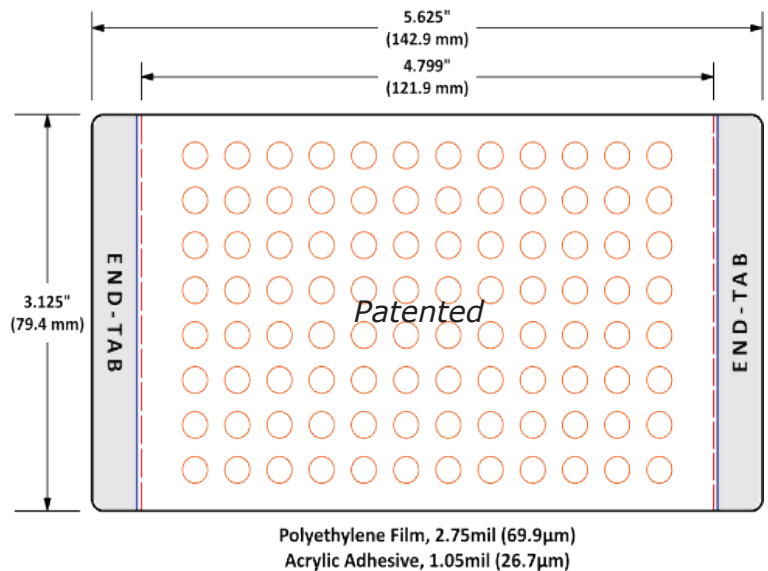

"For a pierceable film with resealability, check out X-Pierce"