

PR1MA™ reagents are available exclusively from MIDSCI™

# PR1MA™ LIFE SCIENCE REAGENTS **CATALOG** POLYMERASES, MASTER MIXES, & DNA LADDERS



## ABOUT MIDSCI™

MIDSCI is pleased to introduce their new **PR1MA** Life Science Reagent Catalog.

Our success has come, by treating our customers and employees with the highest regard. When a company can take care of their team, they will, in turn, take care of their customers the same way.

Our customers are our lifeline, and we will do our very best to bring only quality products to our scientists. In doing so, we will answer our phones with real people who have real answers for you.

- ▶ We ship stocked products in 24 hours.
- ▶ We have exclusive agreements with vendors from all over the world.
- ▶ We will get back to you in 24 hours.
- ▶ We will do all the little things that many have forgotten.

### ALSO AVAILABLE FROM MIDSCI



### PR1MA™ THERMAL CYCLERS



# PR1MA™

## MIDSCI™ Mission

Founded in 1983, MIDSCI has grown to be a prominent name in research supplies, reagents, and equipment thanks to the trust and satisfaction of our customers like you.


As a family-owned certified small business, our Mission is simple:

*Large Enough to Deliver, Small  
Enough to Care.*

We understand the importance of your research and maximizing the effectiveness of your funding dollars. Therefore, we are constantly listening to our customers and working with manufacturers to deliver products that bring not only value, but also process improvement and optimization.

These product advantages are realized in:

- ▶ Sample maximization
- ▶ Reduced waste
- ▶ Time savings
- ▶ Reduced labor costs
- ▶ Or simply easier use by the researcher

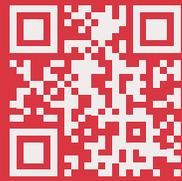

# TABLE OF CONTENTS

<b>1.</b>		<b>ENDPOINT PCR</b>	<b>04</b>
		Taq DNA Polymerase & Master Mix	4
		Hot Start Taq Polymerase & Master Mix	5
		High Fidelity Taq Polymerase & Master Mix	6
		1-Hour Mammalian Genotyping Kit	7
		dNTP Sets & Mixes	7
		One-Step RT-PCR Kit	8
<b>2.</b>		<b>REAL TIME PCR</b>	<b>09</b>
		qMAX™ Green qPCR Mix	9
		qMAX™ Probe qPCR Mix	9
		qMAX™ Probe One-Step RT-qPCR Kit	10
		qMAX™ Probe One-Step Multiplex RT-qPCR Kit	10
		qMAX™ Green One-Step RT-qPCR Kit	11
<b>3.</b>		<b>cDNA SYNTHESIS</b>	<b>12</b>
		qMAX™ cDNA Synthesis Kit	12
		qMAX™ First Strand cDNA Synthesis Kit	12
<b>4.</b>		<b>ELECTROPHORESIS</b>	<b>13</b>
		SmartCheck™ DNA Ladders	13
		SmartGlow™ DNA Stains	13
		SmartCheck™ Protein Markers	14

## Taq DNA Polymerase

- Two-tube formulation: Taq DNA Polymerase & 5X buffer.
- Improved template affinity and solubility for higher enzyme activity and greater yields.
- Exhibits 5'→3' exonuclease activity.
- High performance with gDNA & GC- & AT-rich sequences.

### Each 500 unit kit includes:

- 1 tube of Taq DNA Polymerase (5 U/μL)
- 4 tubes of 5X buffer (PCR enhancers & 15 mM MgCl<sub>2</sub>)

ITEM #	DESCRIPTION
PR1000-500	Taq DNA Polymerase, 500 units
PR1000-1000	Taq DNA Polymerase, 1,000 units
PR1000-6000	Taq DNA Polymerase, 6,000 units

\*dNTPs not included but available separately (see page 7)\*



PR1000-500

### Applications:

- Colony PCR
- Fast PCR
- Genotyping
- High-throughput PCR
- Routine PCR

## Taq 2X Master Mix

- Single-tube formulation: Ready to use 2X Master Mix.
- High performance with gDNA & GC-rich sequences.

### Each mix includes:

- Taq DNA Polymerase
- Buffer [PCR Enhancers, 2 mM dNTPs, 6 mM MgCl<sub>2</sub>, inert red gel loading dye for direct gel loading (optional)]

### Applications:

- Colony PCR
- Fast PCR
- Genotyping
- High-throughput PCR
- Routine PCR



PR1001-R-200

ITEM #	DESCRIPTION
PR1001-200	Taq 2X Master Mix, 200 reactions (4 x 1.25 mL)
PR1001-1000	Taq 2X Master Mix, 1,000 reactions (20 x 1.25 mL)
PR1001-R-200	Taq 2X Red Dye Master Mix, 200 reactions (4 x 1.25 mL)
PR1001-R-1000	Taq 2X Red Dye Master Mix, 1,000 reactions (20 x 1.25 mL)

## ■ Hot Start Taq DNA Polymerase

- Two-tube formulation: Hot Start Taq DNA Polymerase & 5X buffer
- Allows for reaction set-up at room temperature without the risk of non-specific amplification
- Optimized buffer for fast cycling & reproducibility
- High performance with gDNA & GC-rich sequences

### Each 500 unit kit includes:

- 1 tube of Hot Start Taq DNA Polymerase (5 U/μL)
- 4 tubes of 5X buffer (PCR enhancers & 15 mM MgCl<sub>2</sub>)



PR1000-HS-200

### Applications:

- Colony PCR
- Fast PCR
- GC-rich DNA Amplification
- Low Copy PCR
- Multiplex PCR

ITEM #	DESCRIPTION
PR1000-HS-500	Hot Start Taq DNA Polymerase, 500 units
PR1000-HS-1000	Hot Start Taq DNA Polymerase, 1,000 units
PR1000-HS-6000	Hot Start Taq DNA Polymerase, 6,000 units

\*dNTPs not included but available separately (see page 7)\*

## ■ Hot Start Taq 2X Master Mix

- Single-tube formulation: Ready to use 2X Master Mix.
- Allows for reaction set-up at room temperature without the risk of non-specific amplification
- High performance with gDNA & GC-rich sequences

### Each mix includes:

- Taq DNA Polymerase
- Buffer (PCR Enhancers, 2 mM dNTPs, 6 mM MgCl<sub>2</sub>, inert red gel loading dye for direct gel loading (optional))

### Applications:

- Colony PCR
- Fast PCR
- GC-rich DNA Amplification
- Low Copy PCR
- Multiplex PCR



PR1001-HSR-200

ITEM #	DESCRIPTION
PR1001-HS-200	Hot Start Taq 2X Master Mix, 200 reactions (4 x 1.25 mL)
PR1001-HS-1000	Hot Start Taq 2X Master Mix, 1,000 reactions (20 x 1.25 mL)
PR1001-HSR-200	Hot Start Taq 2X Red Dye Master Mix, 200 reactions (4 x 1.25 mL)
PR1001-HSR-1000	Hot Start Taq 2X Red Dye Master Mix, 200 reactions (20 x 1.25 mL)



## High Fidelity Polymerase

- Two-tube formulation: High Fidelity DNA Polymerase & 5X buffer
- 50X the fidelity of wild-type Taq Polymerase
- Exhibits 3'→5' proofreading activity
- High performance with gDNA & GC- & AT-rich sequences
- Ideal for cloning (blunt end), mutagenesis, & microarray applications

### Each 200 unit kit includes:

- 1 tube of High Fidelity Polymerase (2 U/μL)
- 3 tubes of 5X buffer (PCR enhancers, 5 mM dNTPs, & 15 mM MgCl<sub>2</sub>)

### Applications:

- Blunt End Cloning
- Colony PCR
- Fast PCR
- GC- & AT-Rich DNA Amplification
- Low Copy PCR

ITEM #	DESCRIPTION
PR1000-HF-200	High Fidelity DNA Polymerase, 200 units
PR1000-HF-1000	High Fidelity DNA Polymerase, 1,000 units

## High Fidelity 2X Master Mix

- Single-tube formulation: Ready to use 2X Master Mix.
- Ideal for blunt end cloning
- 100X the fidelity of wild-type Taq
- High Fidelity Hot-Start Master Mix available for TA cloning (10X fidelity of wild-type Taq)
- Exhibits 3'→5' proofreading activity

### Each mix includes:


- High Fidelity Polymerase
- Buffer [PCR Enhancers, 2 mM dNTPs, 6 mM MgCl<sub>2</sub>, hot-start antibodies(optional)]

### Applications:


- Blunt End / TA Cloning
- Colony PCR
- Fast PCR
- GC- & AT-Rich DNA Amplification
- Low Copy PCR



PR1001-HFHS-200

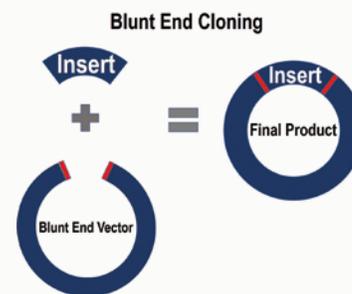

Blunt End Cloning  
(PR1000-HF, PR1001-HF)

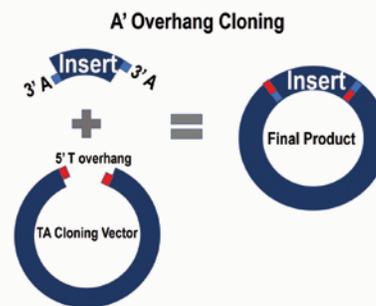

TA Cloning  
(PR1001-HFHS)



ITEM #	DESCRIPTION
PR1001-HF-200	High Fidelity Master Mix, 200 reactions (2 x 1.25 mL)
PR1001-HF-500	High Fidelity Master Mix, 500 reactions (5 x 1.25 mL)
PR1001-HFHS-200	High Fidelity Hot Start Master Mix, 200 reactions (2 x 1.25 mL)
PR1001-HFHS-500	High Fidelity Hot Start Master Mix, 500 reactions (5 x 1.25 mL)

## 1 Hour Mammalian Genotyping Kit

- Three-tube formulation: Lysis buffer, Protease buffer, & Hot Start Taq Red Dye Master Mix
- DNA extraction & amplification in as little as 1 hour
- Designed for use with tail snips, ear punches, blood, FFPE tissue sections, buccal swabs, hair, & solid animal tissues
- Included Hot Start Taq Red Dye Master Mix for reaction set-up at room temperature & direct gel loading

### Each 80 reaction kit includes:

- 1 tube of 5X Lysis Buffer
- 1 tube of 10X Protease Buffer
- 2 tubes of Hot Start Taq 2X Red Dye Master Mix

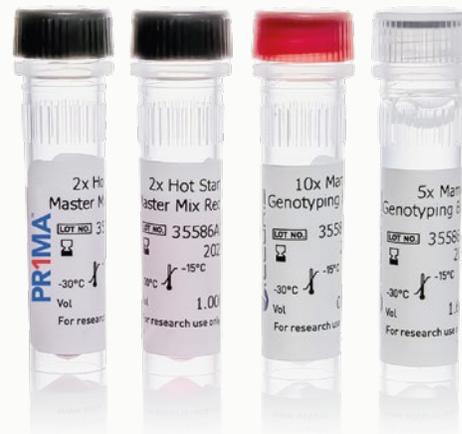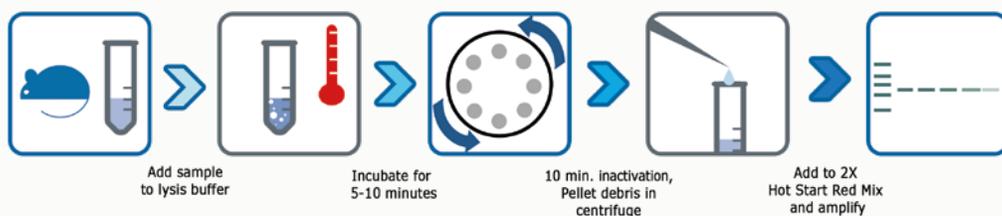

PR1300-MG-80

### Applications:

- DNA Extraction & Amplification from mammalian tissue, blood buccal swabs, & hair

#### ITEM #

#### DESCRIPTION

PR1300-MG-80	1 Hour Mammalian Genotyping Kit, 80 reactions
PR1300-MG-400	1 Hour Mammalian Genotyping Kit, 400 reactions

## dNTP Sets & Mixes

- Supplied as a ready-to-use 40 mM mix or as a set of 4 individual solutions of 100 mM.
- Free of impurities and inhibitors that reduce sensitivity and yield
- Free of DNase, RNase, Nickase, & Protease activity (>99% pure, purified by HPLC)
- 24 month shelf life

### Each 0.5 mL mix includes:

- 1 tube of 40 mM dNTP Mix

### Each 4 x 250 µL set includes:

- 1 tube of 100 mM dATP
- 1 tube of 100 mM dCTP
- 1 tube of 100 mM dGTP
- 1 tube of 100 mM dTTP



PR3101-S-1

#### ITEM #

#### DESCRIPTION

PR3040-M-1	40 mM dNTP Mix (10 mM/20 µmol each), 0.5 mL
PR3040-M-4	40 mM dNTP Mix (10 mM/80 µmol each), 2 mL
PR3101-S-1	100 mM dNTP Set, dATP, dGTP, dTTP, dCTP, 250 µL /25 µmol of each dNTP
PR3101-S-4	100 mM dNTP Set, dATP, dGTP, dTTP, dCTP, 1 mL/100 µmol of each dNTP

### Applications:

- End-Point PCR
- Real-Time PCR
- cDNA Synthesis


- Starting Material: **RNA**
- Two-tube formulation: 2X One-Step Mix & 20X RTase Blend
- Convenient cDNA synthesis and subsequent PCR in a single tube
- Sensitive: Optimized chemistry for detection of low copy transcripts
- Ideal for multiplex PCR & gene expression analysis
- Hot Start Taq Polymerase included to prevent nonspecific amplification and allow sample preparation at room temperature

**Each 50 reaction kit includes:**

- 1 tube of 20X thermostable RTase (RNase inhibitor included)
- 1 tube of 2X One-Step Mix (PCR enhancers, 2 mM dNTPs, & 6 mM MgCl<sub>2</sub>)



PR1100-50

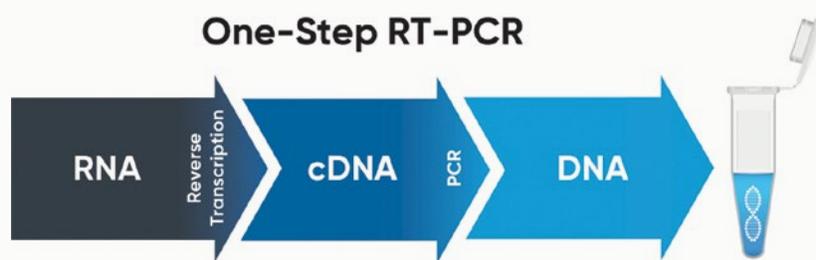

**Applications:**

- Multiplex RT-PCR
- Gene Expression Analysis

ITEM #	DESCRIPTION
PR1100-50	One-Step RT-PCR Kit, 50 reactions
PR1100-100	One-Step RT-PCR Kit, 100 reactions

Copyright © 2024 Benchmark Scientific





## qMAX™ Green qPCR Mix

- Single-tube formulation: Ready to use 2X Master Mix (Includes Hot Start Taq DNA Polymerase)
- Proprietary intercalating green fluorescent dye for minimal PCR inhibition & optimal fluorescence
- Inert blue tracking dye to assist in visualization during pipetting
- Early Ct values and detection across a broad dynamic range

### Each mix includes:

- PR1MA Hot Start Taq DNA Polymerase
- qMAX Green (proprietary fluorescent dye)
- Inert blue tracking dye
- dNTPs
- Optimized buffer for maximum efficiency & sensitivity



### Applications:

- Gene Expression Analysis
- Fluorescent Detection of DNA/cDNA



PR2000-N-100

ITEM #	DESCRIPTION
PR2000-N*-100	qMAX™ Green qPCR mix with blue tracking dye, No ROX, 100 reactions (1 x 1.0 mL)
PR2000-N*-500	qMAX™ Green qPCR mix with blue tracking dye, No ROX, 500 reactions (5 x 1.0 mL)
PR2000-N*-1000	qMAX™ Green qPCR mix with blue tracking dye, No ROX, 1,000 reactions (10 x 1.0 mL)

\* To specify ROX levels when ordering: Use -N for No ROX, -L for Low ROX, & -H for High ROX  
Bulk packaging available upon request



## qMAX™ Probe qPCR Mix


- Single-tube formulation: Ready to use 2X Master Mix (Includes Hot Start Taq DNA Polymerase)
- Compatible with popular hydrolysis & beacon probes
- Ideal for cDNA, gDNA, & low-copy number templates
- Compatible with standard/fast cycling protocols.

### Each mix includes:

- PR1MA Hot Start Taq DNA Polymerase
- dNTPs
- Optimized buffer for maximum efficiency & sensitivity



PR2001-N-100

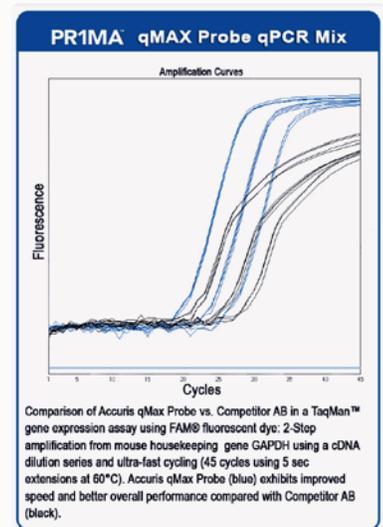

### Applications:

- Multiplex RT-qPCR
- Probe-based detection of DNA/cDNA
- Gene Expression Analysis

ITEM #	DESCRIPTION
PR2001-N*-100	qMAX™ Probe qPCR Mix, No ROX, 100 reactions (1 x 1.0 mL)
PR2001-N*-500	qMAX™ Probe qPCR Mix, No ROX, 500 reactions (5 x 1.0 mL)
PR2001-N*-1000	qMAX™ Probe qPCR Mix, No ROX, 1,000 reactions (10 x 1.0 mL)

\* To specify ROX levels when ordering: Use -N for No ROX, -L for Low ROX, & -H for High ROX  
Bulk packaging available upon request



### ■ Probe One-Step RT-qPCR Kit

- Starting Material: **RNA**
- Two-tube formulation: 2X Probe One-Step Mix & 20X RTase Blend (Includes Hot Start Taq Polymerase)
- Convenient cDNA synthesis and subsequent qPCR in a single tube
- Ideal for gene expression analysis & **single** pathogen detection

#### Each 100 reaction kit includes:


- 1 tube of 20X thermostable MMLV derived RTase (RNase inhibitor included)
- 1 tube of 2X Probe One-Step Mix
- Includes PR1MA Hot Start Taq DNA Polymerase
- dNTPs
- Optimized buffer for maximum efficiency & sensitivity

#### Applications:

- Single Pathogen Detection
- Gene Expression Analysis



PR2121-N-100

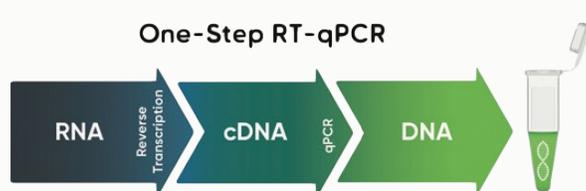

### ■ Probe One-Step Multiplex RT-qPCR Kit

- Starting Material: **RNA**
- Two-tube formulation: 2X Probe One-Step Mix & 20X RTase Blend (Includes Hot Start Taq Polymerase)
- Convenient cDNA synthesis and subsequent qPCR in a single tube
- Ideal for gene expression analysis & **multiple** pathogen detection

#### Each 100 reaction kit includes:


- 1 tube of 20X thermostable MMLV derived RTase (RNase inhibitor included)
- 1 tube of 2X Probe One-Step Mix
- Includes PR1MA Hot Start Taq DNA Polymerase
- dNTPs
- Optimized buffer for maximum efficiency & sensitivity

#### Applications:

- Multiplex RT-qPCR
- Gene Expression Analysis



PR2122-N-100

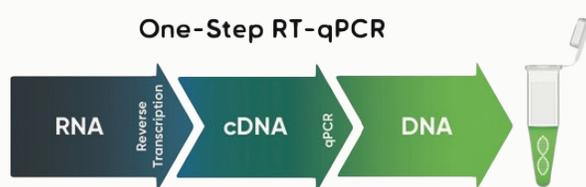

ITEM #	DESCRIPTION
PR2121-N*-100-500-1000	qMAX™ Probe One-Step RT-qPCR Kit, No ROX 100/500/1,000 reactions
PR2122-N*-100-500-1000	qMAX™ Probe One-Step Multiplex RT-qPCR Kit, No ROX 100/500/1,000 reactions

\* To specify ROX levels when ordering: Use -N for No ROX, -L for Low ROX, & -H for High ROX  
Bulk packaging available upon request



- Starting Material: **RNA**
- Two-tube formulation: 2X Green One-Step Mix & 20X RTase Blend (Includes Hot Start Taq Polymerase)
- Proprietary intercalating green fluorescent dye for minimal PCR inhibition & optimal fluorescence
- Earlier Ct values and detection across a broad dynamic range
- Ideal for cDNA, gDNA, & low-copy number templates




PR2120-N-100

### Each 100 reaction kit includes:

- 1 tube of 20X thermostable MMLV derived RTase (RNase inhibitor included)
- 1 tube of 2X Green One-Step Mix
- Includes PR1MA Hot Start Taq DNA Polymerase
- qMAX Green (proprietary fluorescent dye)
- dNTPs
- Optimized buffer for maximum efficiency & sensitivity

### Applications:

- Gene Expression Analysis
- Fluorescent Detection of DNA/cDNA

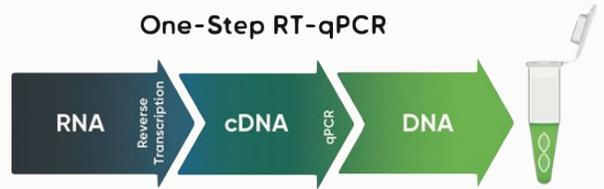

REAL TIME PCR

ITEM #	DESCRIPTION
PR2120-N*-100	qMAX™ Green One-Step RT-qPCR Kit, 100 reactions
PR2120-N*-500	qMAX™ Green One-Step RT-qPCR Kit, 500 reactions
PR2120-N*-1000	qMAX™ Green One-Step RT-qPCR Kit, 1,000 reactions

\* To specify ROX levels when ordering: Use -N for No ROX, -L for Low ROX, & -H for High ROX  
Bulk packaging available upon request



## ■ cDNA Synthesis Kit

- Starting Material: **RNA**
- Two-tube formulation: 20X RTase & 5X Buffer
- Ideal for 4 pg - 0.5 µg of input RNA

### Each 25 reaction kit includes:


- 1 tube of 20X RTase (10 U/µL; potent RNase inhibitor included)
- 1 tube of 5X buffer (oligo dT primers & random hexamers included)

### Applications:

- cDNA Library Generation
- Microarrays
- Gene Expression Analysis
- 2-Step RT-PCR/RT-qPCR



PR2100-C-25

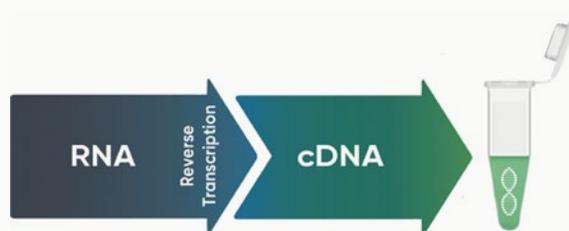

## ■ First Strand cDNA Synthesis Flex Kit

- Starting Material: **RNA**
- Four-tube formulation for greater flexibility in assay design
- Ideal for 10 pg - 2.0 µg of input RNA

### Each 50 reaction kit includes:


- 1 tube of RTase (200 U/µL; Potent RNase inhibitor included)
- 1 tube of 5X buffer (5 mM dNTPs, 15 mM MgCl<sub>2</sub>, stabilizers, & enhancers)
- 1 tube of random hexaprimers
- 1 tube of oligo dT primers

### Applications:

- cDNA Library Generation
- Microarrays
- Gene Expression Analysis
- 2-Step RT-PCR/RT-qPCR



PR2110-50

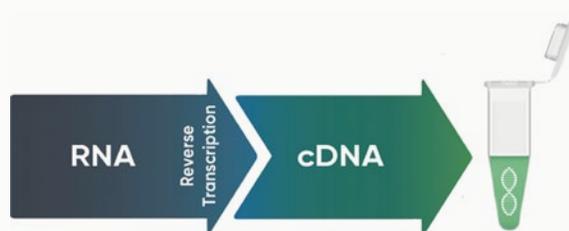

### ITEM #

### DESCRIPTION

PR2100-C-25-100-250

qMAX™ cDNA Synthesis Kit, 25-100-250 reactions

PR2110-50-200

qMAX™ First Strand cDNA Synthesis Flex Kit, 50/200 reactions




## SmartCheck™ DNA Ladders


- Ready to use formulation: includes loading buffer & tracking dye
- Three ranges: 50 bp, 100 bp, & 1kb
- Purified by multistep chromatography to produce clear, distinct bands with increased intensity for easy size determination

### Each 100 lane kit includes:


- 1 x 500 µL tube of SmartCheck DNA Ladder (0.1 mg/mL) supplied in buffer

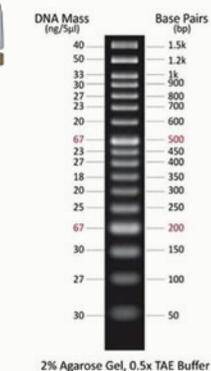

PR4005-100

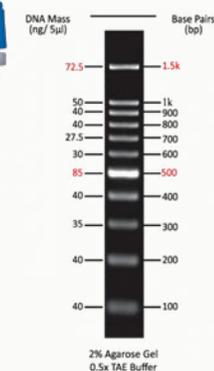

PR4010-100

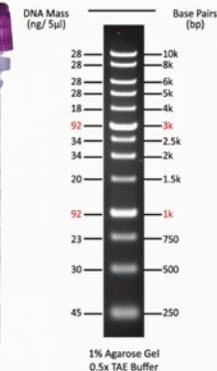

PR4100-100

ITEM #	DESCRIPTION
PR4005-100-500	SmartCheck™ 50bp DNA Ladder, 100-500 lanes
PR4010-100-500	SmartCheck™ 100bp DNA Ladder, 100-500 lanes
PR4100-100-500	SmartCheck™ 1kb DNA Ladder, 100-500 lanes

ELECTROPHORESIS



## SmartGlow™ DNA Stains

- Non-mutagenic, environmentally safe
- Excitation by UV (290nm) or blue light (490nm)
- Two types available: PS (Pre-Stain) & LD (Loading Dye)
- Emits green fluorescence when bound to ssDNA, dsDNA, & RNA

### Each Pre-Stain tube:

- 1 tube of 20,000X DNA Pre-Stain

### Each Loading Dye tube:


- 1 tube of 6X DNA Loading Dye



E4500-PS



E4500-LD

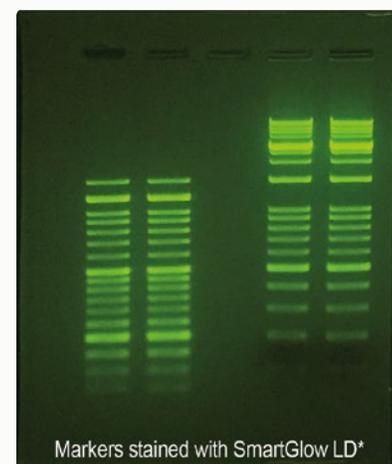

Markers stained with SmartGlow LD\*



ITEM #	DESCRIPTION
BE4500-PS	SmartGlow™ Pre-Stain 1.0 mL (20,000X)
BE4500-LD	SmartGlow™ Loading Dye, 1.0 mL (6X)



- Ready to use formulation
- 13 pre-stained protein standards
- Molecular weight range: 5 - 245 kDa
- Recommended load volume: 5  $\mu$ L / well

### Each kit includes:

- 2 x 250  $\mu$ L of Tri-Color Protein Marker (0.1-0.4 mg/mL) in buffer (20 mM Tris-Phosphate, 2% SDS, 0.2 mM Dithiothreitol, 3.6 M Urea, & 15% Glycerol)

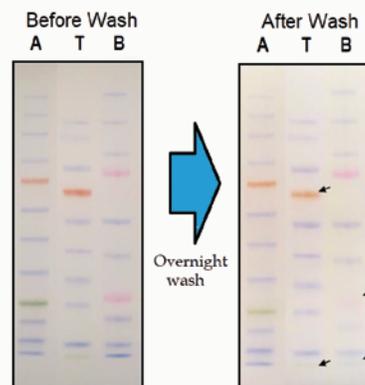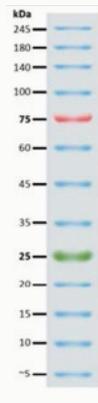

Customer testing shows that Accuris SmartCheck™ Protein Markers (A), perform better than competitor brands (T and B) with less fading of bands after overnight washing. Bands were transferred to nitrocellulose membrane and washed with TBS containing 0.2% Tween 20 overnight.

ITEM #	DESCRIPTION
PR5013-TC	SmartCheck™ Tri-color Protein Marker, 13 bands 5-245 kDa

Copyright © 2024 MIDSCI™



CONTINUOUS IMPROVEMENT

ALSO AVAILABLE FROM MIDSCI™

## AutoMATE™ 96

Microplate Pipetting Workstation



## isoPURE™

Magnetic Purification System




## MYVOLT™ TOUCH



## SmartDrop™



## MY GEL™ InstaView

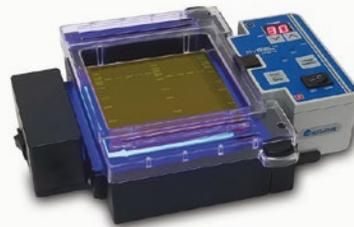

**MIDSCI™**  
Laboratory Equipment & Supplies

*Large Enough to Deliver, Small Enough to Care.*



**CONTACT US**

115 Cassens Court  
Fenton, MO, 63026  
United States

800-227-9997  
Customer Service: [custserv@midsci.com](mailto:custserv@midsci.com)  
[www.midsci.com](http://www.midsci.com)

[www.midsci.com](http://www.midsci.com)